


# Danis Procedure Pack

SX-ELLA Stent Danis

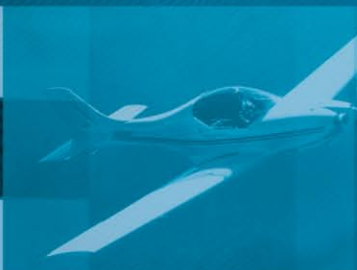

Removable stent system for treatment of bleeding esophageal varices

Easy implantation without the need for endoscopy / X-Ray control

Secure hemostasis through direct compression of varices

Oral intake maintained from time of implantation

## Basic description

### Procedure Pack with SX-ELLA Stent Danis

Self-expandable nitinol stent preloaded in a ready-to-use delivery system  
Complete procedure pack including all items required for implantation

#### Indications:

Acute bleeding from esophageal varices

#### Contraindications:


Esophageal strictures  
Recent radiotherapy of the esophagus  
Tumors of the upper respiratory tract, esophagus or stomach  
Upper GIT bleeding from other causes (e.g. perforation)

#### Features / Benefits:


- ✓ **Variable pitches in stent braiding** conform to esophageal peristalsis, reducing the risk of stent migration
- ✓ **Standardized compression of the varices** results in effective hemostasis
- ✓ **Atraumatic edges**
- ✓ **Radiopaque markers** at both stent ends and at the mid-point
- ✓ **Retrieval loops** with gold markers at both stent ends
- ✓ **Easy endoscopic stent removal** after 7 days

#### Benefits of the method


- ✓ Readily implantable without endoscopy / X-Ray control, also in emergency situations
- ✓ Patient can continue an oral diet after stent implantation

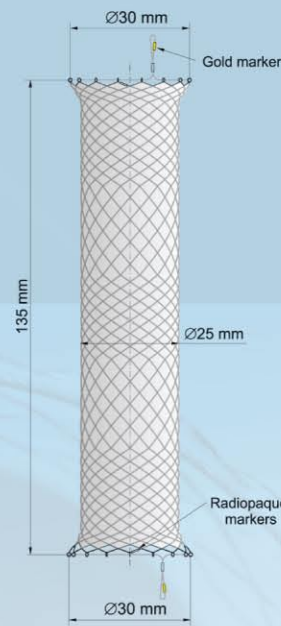

#### Delivery system

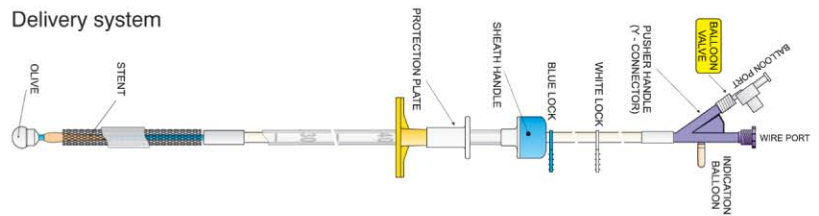

#### Delivery system - Gastric balloon deployed

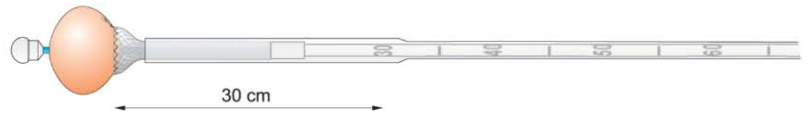

#### Additional description:

The procedure pack consists of a plastic carry case containing all items required for stent implantation:  
Delivery system pre-loaded with stent, guide wire, plastic syringes, mouthpiece, kidney dish, slobber cloth, waste bag with clip, sterile gloves, instruction for use, patient's card and self-adhesive product label.

Nominal (relaxed) diameter of stent body is 25 mm, nominal diameter of throats is 30 mm, and nominal stent length is 135 mm.

The delivery system 8.2 mm (24.5 F) / 6.2 mm (19 F) is available in a standard length of 64 cm.

REF	SX-ELLA Stent Danis		Delivery system	
	Nominal diameter [mm]	Nominal length [mm]	Active length [cm]	Outer diameter [F]
019-08-25-135	25	135	64	24.5 F / 19 F

#### Manufacturer:

**ELLA-CS, s.r.o.**

Milady Horakove 504, 500 06 Hradec Kralove  
Czech Republic

Tel: +420 495 279 111 Fax: +420 495 265 655  
E-mail: info@ellacs.eu URL: www.ellacs.eu