



Enterella

SX-ELLA Stent Colorectal

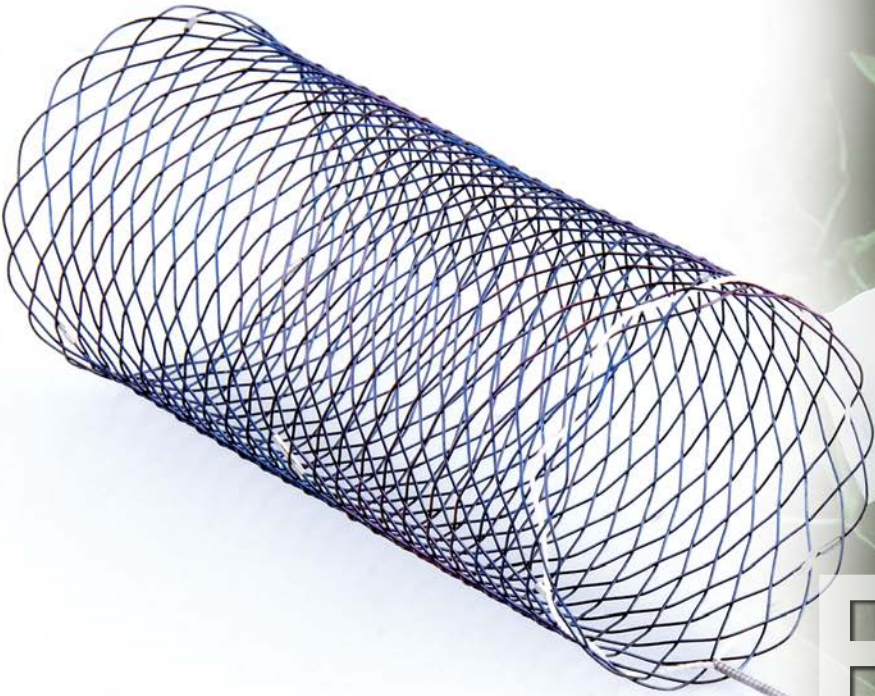

Enterella

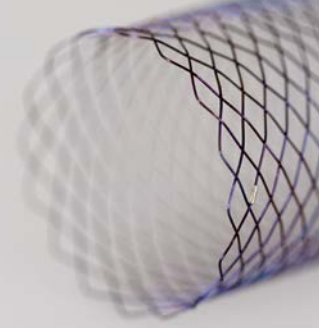

Easy deployment

Excellent flexibility

Atraumatic stent ends

Smooth alignment around flexures



Basic description

SX-ELLA Stent Colorectal – Enterella

Stent system for insertion under fluoroscopic or endoscopic guidance
Self-expandable nitinol stent preloaded in a pull-back delivery system

Indications

Temporary or definitive palliative treatment of malignant strictures affecting the rectum, sigmoid colon or distal part of descending colon
Treatment of coloenteral and / or colovesical fistulas (covered stent)

Contraindications

- Bowel perforation
- Peritoneal abscess
- Severe ascites
- Severe coagulopathy
- Enteral ischemia
- Placement of intestinal obstruction / stricture outside the anatomic structures defined in "Indications"

Features / Benefits


- ☑ Atraumatic stent ends / Effective prevention of trauma to the bowel wall
- ☑ Excellent flexibility / Precise alignment around anatomical curves
- ☑ Radiopaque markers at both stent ends and at the midpoint / Improved visibility; accurate stent positioning
- ☑ Plastic repositioning loop at aboral stent end / Possibility of stent repositioning within five days after implantation
- ☑ Optional stent covering / Prevention of tumor in-growth; occlusion of fistulas

Additional description


15 F (5.08 mm) delivery system available in a standard length of 95 cm (135 cm delivery system is available on request, covered stent requires 18 F (6.0 mm) delivery system).

Radiopaque, atraumatic olive tip of the delivery system.
Recommended guide wire: 0.035" (0.89 mm) ultra stiff.



Excellent flexibility

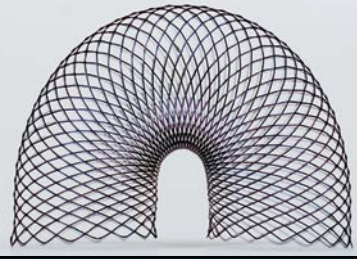

Plastic repositioning loop

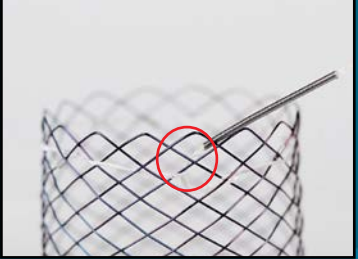

Atraumatic stent ends

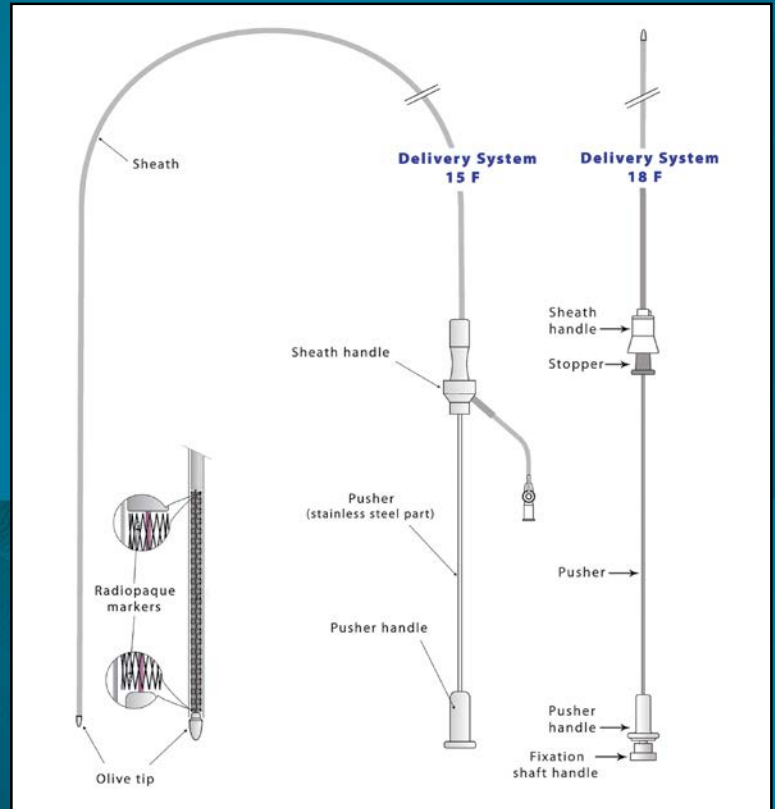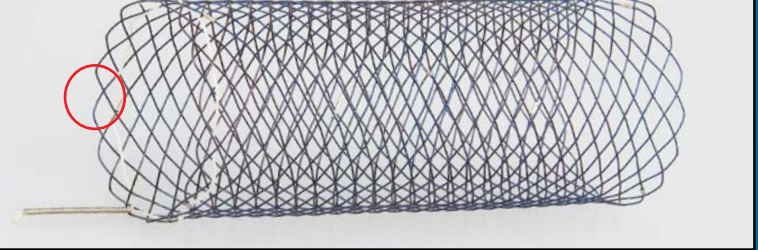

Available sizes

REF No. Uncovered version	REF No. Covered version	Nominal (unconstrained) diameter [mm]	Nominal (unconstrained) length [mm]
019-02-22-082	019-02-22-082-C18	22	85
019-05-22-090	019-05-22-090-C18	22	90
019-05-22-113	019-05-22-113-C18	22	113
019-05-22-135	019-05-22-135-C18	22	135
019-05-25-082	019-05-25-082-C18	25	82
019-05-25-090	019-05-25-090-C18	25	90
019-05-25-113	019-05-25-113-C18	25	113
019-05-25-135	019-05-25-135-C18	25	135
019-05-30-075	019-05-30-075-C18	30	75
019-05-30-088	019-05-30-088-C18	30	88
019-05-30-112	019-05-30-112-C18	30	112
019-05-30-123	019-05-30-123-C18	30	123
019-05-30-136	019-05-30-136-C18	30	136

Easy deployment



Distributor:

Manufacturer:

ELLA-CS, s.r.o.

Milady Horakove 504/45, Trebes,
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655
E-mail: info@ellacs.eu URL: www.ellacs.eu