


Flexella Plus

SX-ELLA Stent Esophageal



Flexella Plus

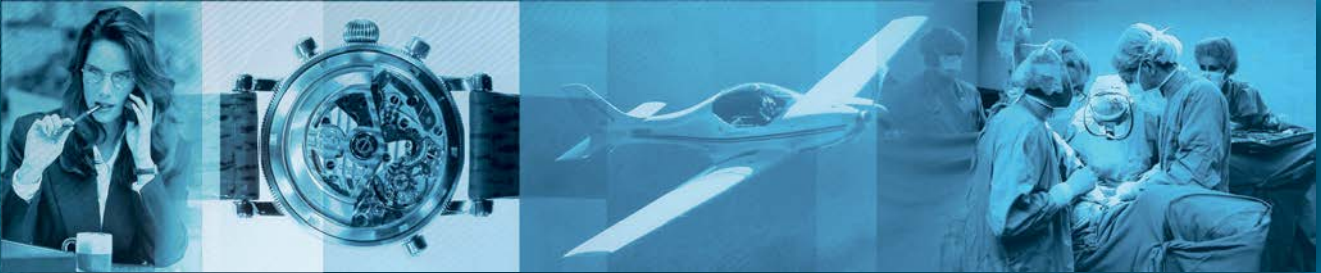


Durable silicone covering

Removable esophageal stent

Anti-migration and atraumatic flares

Excellent flexibility

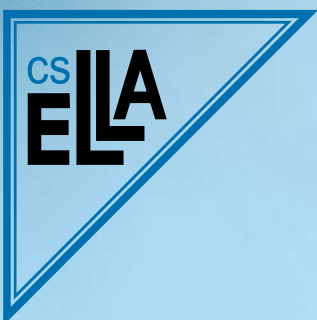

Basic description

SX-ELLA Stent Esophageal – Flexella Plus

Flexella Plus is designed for palliation of malignant esophageal strictures, occlusion of esophagorespiratory fistulas and benign strictures resistant to standard therapy.

Stent system for fluoroscopic and / or endoscopic insertion
Self-expandable fully covered nitinol stent preloaded in PULL delivery system

Features / Benefits





- ☑ 2 types of delivery system: PULL and PUSH
- ☑ PUSH delivery system / Strictures located in upper part of esophagus
- ☑ PULL delivery system / Splittable olive tip / Automatic detachment while stent starts releasing / Prevention of olive entrapment / Safe removal of delivery system
- ☑ Retrieval loop made of medical grade alloy at proximal stent end / High mechanical durability and acid-resistance / Centripetal constriction of stent flare / Suitable for extraction by means of ELLA Extractor
- ☑ Radiopaque markers / Accurate stent positioning
- ☑ Anti-migration design / Flared design / Reduced migration rates
- ☑ Atraumatic stent ends / Reduction of esophageal wall irritation and hyperplasia
- ☑ Excellent flexibility / Smooth alignment and adequate stent expansion along anatomical curves
- ☑ Durable silicone stent covering / Prevention of tumor in-growth, occlusion of esophagorespiratory fistulas

Additional description

PUSH delivery system (model 06S) 4.7 mm / 9.4 mm (14 F / 28 F), available in standard length of 55 cm

PULL delivery system (model 07S) 6 mm / 9.4 mm (18 F / 28 F), available in standard length of 75 cm

Recommended guide wire: 0.035" (0.89 mm) ultra stiff / Minimum length 220 cm

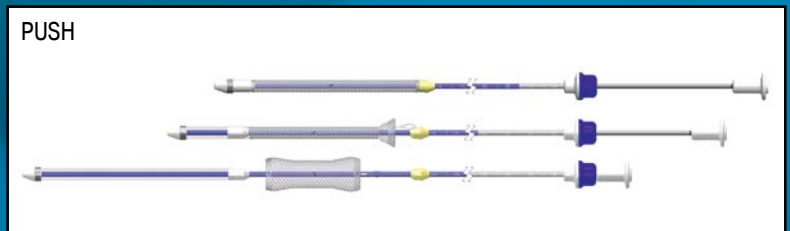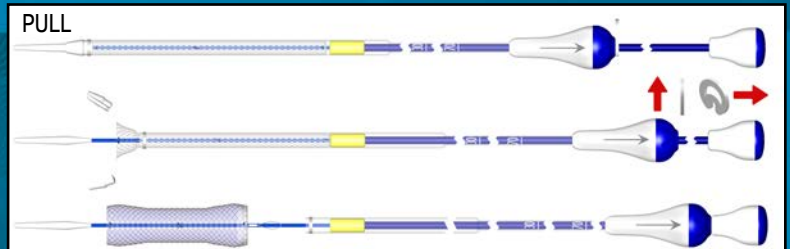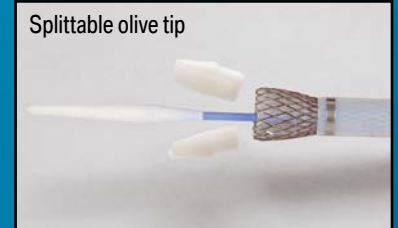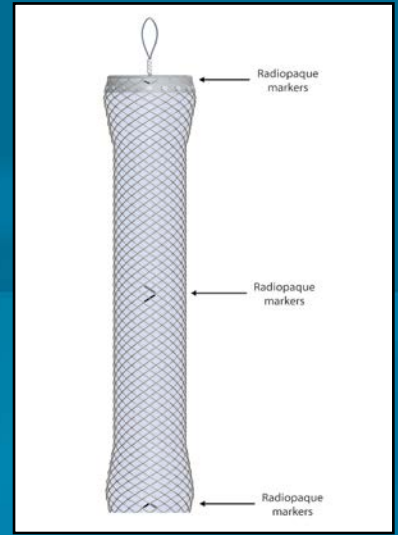

Available sizes

REF No.	Stent body diameter [mm]	Stent flares diameter [mm]	Stent nominal length [mm]	Anti-reflux valve
019-06S-20-085	20	25	85	No
019-06S-20-110	20	25	110	No
019-07S-20-085	20	25	85	No
019-07S-20-110	20	25	110	No
019-07S-20-135	20	25	135	No
019-07S-20-150	20	25	150	No
019-07S-20-085-0	20	25	85	Yes
019-07S-20-110-0	20	25	110	Yes
019-07S-20-135-0	20	25	135	Yes
019-07S-20-150-0	20	25	150	Yes

Distributor:

Manufacturer:

ELLA-CS, s.r.o.

Milady Horakove 504/45, Trebes,
500 06 Hradec Kralove, Czech Republic

Tel.: +420 495 279 111 Fax: +420 495 265 655
E-mail: info@ellacs.eu URL: www.ellacs.eu