


Infinity[®] ERCP sampling device



Looking to collect quality samples during ERCP?

The Infinity® ERCP sampling device is purpose built for collecting substantial and quality samples from strictures in the biliary duct.

The Infinity® device features...

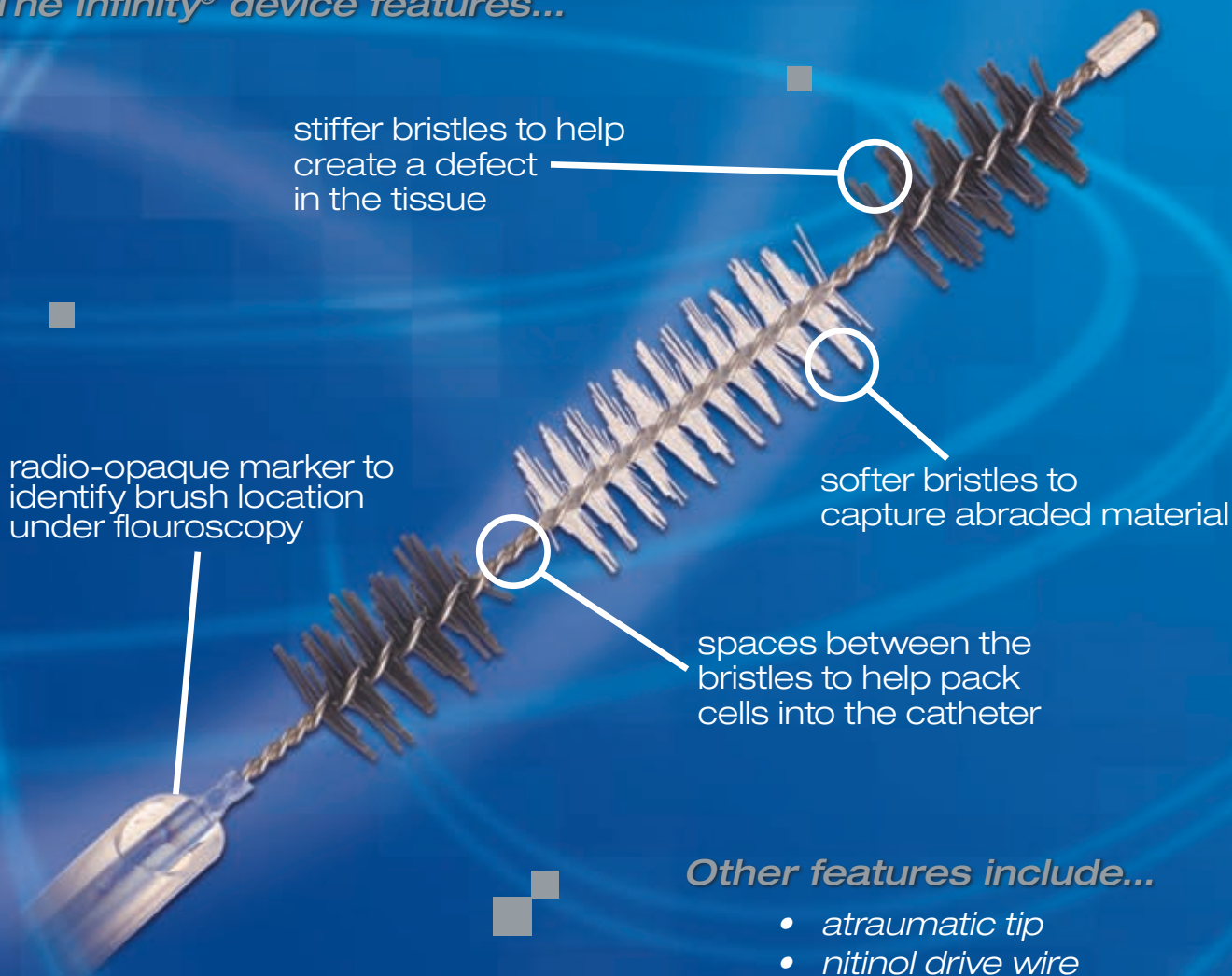

stiffer bristles to help create a defect in the tissue

radio-opaque marker to identify brush location under flouroscopy

softer bristles to capture abraded material

spaces between the bristles to help pack cells into the catheter

Other features include...

- atraumatic tip
- nitinol drive wire
- injection port
- salvage cytology
- long and short guidewire compatibility

For more information on the Infinity® ERCP sampling device, contact US Endoscopy today.

product number	sheath diameter (mm)	length (cm)	sterile	units/box
00711652	3 / 9FR	200	yes	5

Distributed by:



listening...and delivering solutions®

5976 Heisley Road
Mentor, OH 44060 USA

phone +1 440 / 639.4494

fax +1 440 / 639.4495

email global@usendoscopy.com

www.usendoscopy.com

760746 Rev.D