


SmartBand™ Ligation System

Banding is our Business

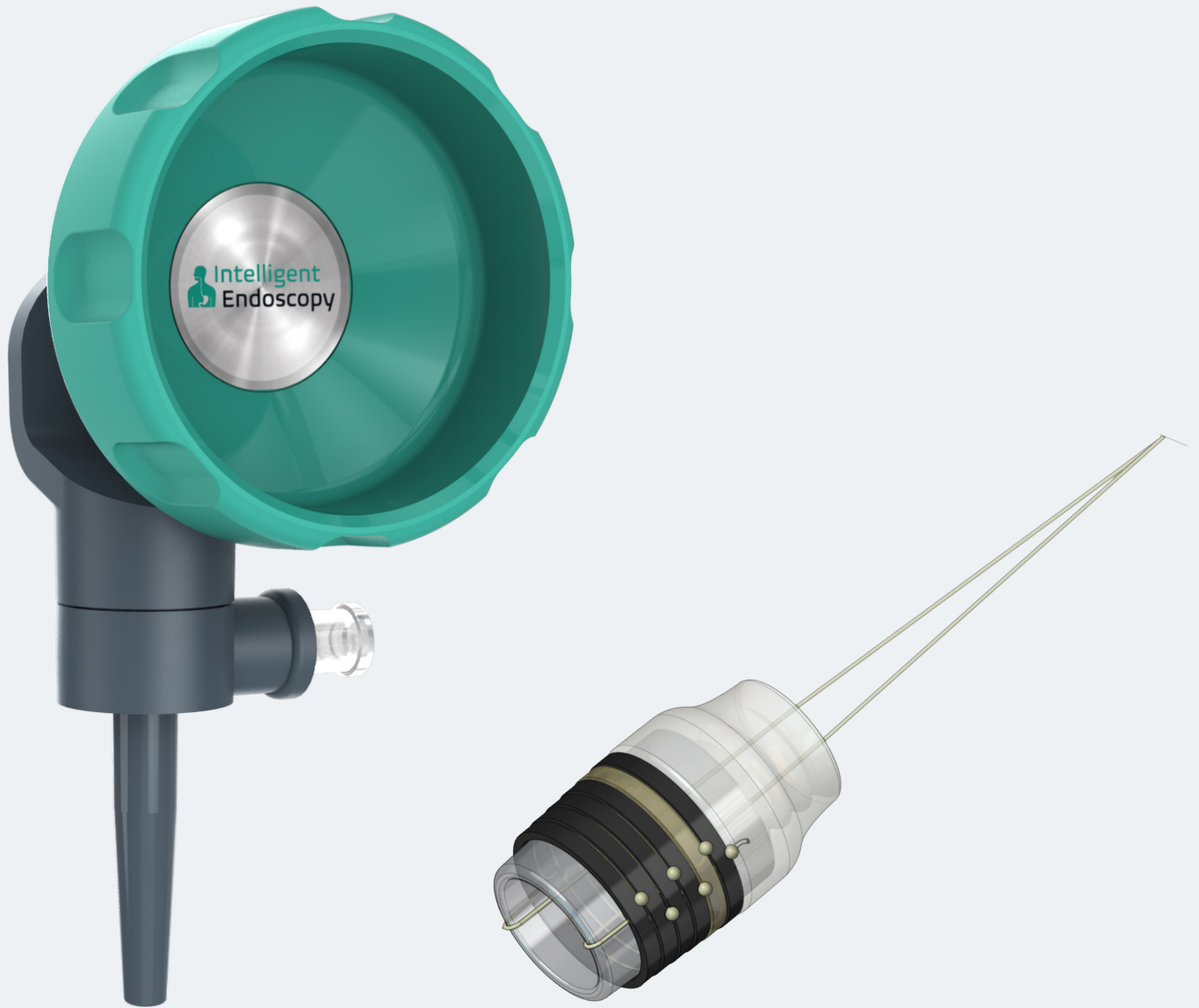

Product Specifications

Used to endoscopically ligate esophageal varices at or above the gastroesophageal junction and to ligate internal hemorrhoids



- High-retention bands contain less than half the latex proteins compared to other banding devices*
- Specialized UV protected packaging maintains band performance for up to 24 months shelf life
- Advanced friction-fit barrel design is compatible with the widest range of upper GI endoscopes from 8.6mm – 11.6mm
- Revolutionary, single-hub ligation handle allows for full-functionality with the simplicity of a one-way rotation for band deployment
- Larger handle stem inner diameter, 9.8Fr, ensures maximum space for suction and flushing capabilities**
- Improve endoscopic view during the procedure by using the dedicated flush port and connector tube for immediate irrigation
- The one-piece loading device utilizes PermaCoil™ technology for reliable device assembly
- Use the standalone SmartBand Ligation Pack when additional bands are required to complete the procedure

SmartBand Ligation Kit

Components: Ligation Handle, Loading Device, Connector Tube and Deployment Cord with Barrel of Bands

Order Number	Reference Part Number	Barrel Diameter (mm)	Number of Bands	Band Material	Shelf Life (months)
IE11002	SLK-6	8.6 - 11.6	6	Latex	24

This device is supplied non-sterile and is disposable, single-use only

SmartBand Ligation Pack

Components: Deployment Cord with Barrel of Bands

Order Number	Reference Part Number	Barrel Diameter (mm)	Number of Bands	Band Material	Shelf Life (months)
IE11001	SLP-6	8.6 - 11.6	6	Latex	24

This device is supplied non-sterile and is disposable, single-use only

*Data on file at Intelligent Endoscopy

**Not compatible with Pentax® endoscopes. Pentax is a registered trademark of Pentax Medical.